



# Field Service Procedure


Part Number: SP00288  
Rev: A  
Date: 19 December 2003  
© 2003 Draeger Medical, Inc.

## AGS Scavenger Installation

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)  
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

## AGS Scavenger Installation

### INSTALLATION PROCEDURE

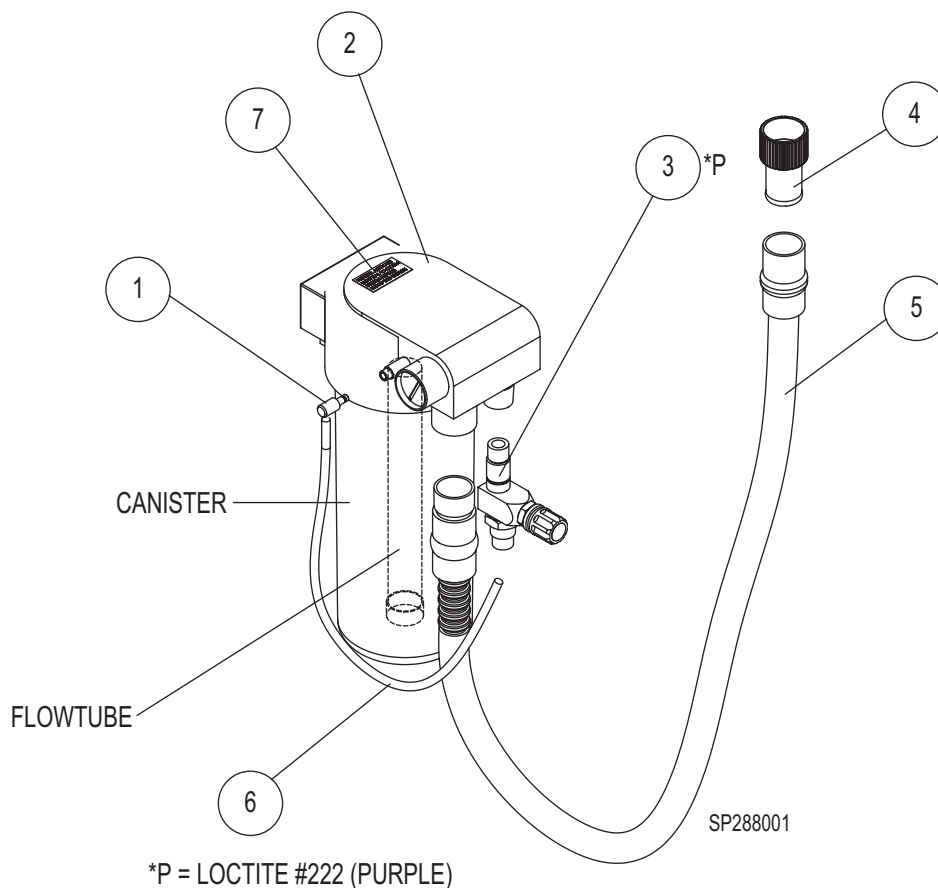

**FIGURE 1. AGS Scavenger Installation**

1. Remove the canister and flow tube from the bottom of the scavenger housing (2).
2. If equipped, apply a small amount of Loctite #222 to the threads of the adjustment valve (3) and screw the valve into the small port of the scavenger (2) until the threads bottom, then back off the valve to the position as shown in Figure 1, but not more than one full turn from bottom position.
3. Position the flow tube so the graduations face toward the front of the machine when the scavenger is mounted.
4. Reinstall the canister.
5. Insert the metal end of the scavenger hose (5) into the large port on the scavenger (2) and twist to secure the hose.

|                                    |
|------------------------------------|
| INSTALLATION PROCEDURE (continued) |
|------------------------------------|

6. Push the 30 mm adapter (4) onto the other end of the scavenger hose (5).
7. Connect the 30 mm adapter to the scavenger port on the underside of the breathing system and twist to secure.
8. If applicable, cut the smaller diameter hose (6) to the desired length for connecting an agent monitor exhaust to the scavenger.
9. Slide the hose onto the barb of the quick disconnect fitting (1) and connect the fitting to the scavenger.
10. Apply label (7) in the approximate position shown in Figure 1.

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)  
[RETURN TO CD-ROM TABLE OF CONTENTS](#)



DrägerService is a division of  
Draeger Medical, Inc.  
3122 Commerce Drive  
Telford, PA 18969  
Tel: (215) 721-5402  
(800) 543-5047  
Fax: (215) 721-5784  
Web: [www.draegermedical.com](http://www.draegermedical.com)  
Printed in the U.S.A.